

IRON (Fe²⁺) IN BHAVANI RIVER WATER: POTENTIOMETRIC ANALYSIS & POLLUTION ASSESSMENT (CPCB STANDARDS)

Exp. No. :


Date :

AIM/ OBJECTIVE:

To estimate the amount of ferrous ion (Fe²⁺) present in the whole of the given solution potentiometrically. A standard solution of potassium dichromate of strengthN is provided.

PRINCIPLE:

Potentiometric titrations depend on measurement of EMF between reference electrode and an indicator electrode (calomel). When a solution of ferrous ion is titrated with a solution of potassium dichromate, the following redox reaction takes place.

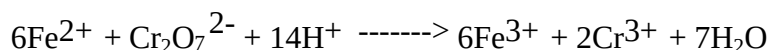

During this titration Fe²⁺ is converted into Fe³⁺ ion and relatively potential of the cell is increases. At the end point, there will be a sharp potential change due to sudden removal of all Fe²⁺ ions. A graph between EMF measured against the volume of K₂Cr₂O₇ added is drawn and the end point is noted from the graph as V₁ ml. The cell is set up by connecting this redox electrode with a calomel electrode as shown below.



MATERIALS REQUIRED:

Potentiometer, Pt electrode, saturated calomel electrode, glass rod, standard K₂Cr₂O₇ solution, ferrous ion (Fe²⁺) solution.


PROCEDURE:

The given ferrous ion solution is made up to 100 ml in a standard flask. 20 ml of this made up solution is pipetted out into a clean 100 ml beaker. About 10 ml of dil. H₂SO₄ and 10 ml of distilled water are added into it. A platinum electrode is dipped into the solution. This electrode is coupled with a saturated calomel electrode and the cell is introduced into potentiometric circuit. The std. K₂Cr₂O₇ solution is taken in the burette and is added in 0.5 ml portions. The emf of the cell is measured after each addition. The addition of K₂Cr₂O₇ is continued even after the end point for further 3 ml. The accurate end point is determined by plotting $\Delta E/\Delta V$ Vs Average volume of K₂Cr₂O₇ added (graph). From the end point, the strength of ferrous ion solution and its amount can be calculated.

OBSERVATIONS/ INFERENCE:**TITRATION OF $K_2Cr_2O_7$ Vs. FERROUS ION SOLUTION****Volume of Ferrous ion solution =ml**

S.No	Volume of $K_2Cr_2O_7$ (ml)	EMF (mV)	ΔE (mV)	ΔV (ml)	$\Delta E / \Delta V$ (mV/ml)	Avg.vol.of $K_2Cr_2O_7$ (ml)
1						
2						
3						
4						
5						
6						
7						
8						
9						
10						
11						
12						
13						
14						
15						
16						
17						
18						
19						
20						
21						
22						
23						
24						
25						
26						
27						

Model graph:

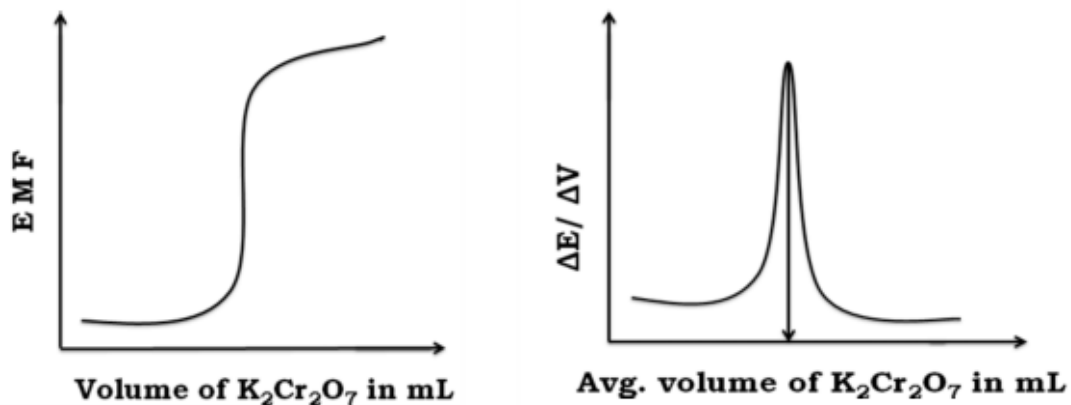

TITRATION OF $K_2Cr_2O_7$ Vs. FERROUS ION SOLUTION

Volume of $K_2Cr_2O_7$ solution V_1 = ml (from graph)

Strength of $K_2Cr_2O_7$ solution N_1 = N

Volume of Ferrous ion solution V_2 = ml

Strength of Ferrous ion solution N_2 = ?

According to the law of volumetric analysis, $V_1N_1 = V_2N_2$

Strength of Ferrous ion (N_2) $N_2 = V_1N_1 / N_2$

=

Strength of Ferrous ion (N_2) = N

Amount of Ferrous ion present in 1 litre of the
given solution

= Eq.Wt. x Stength of Ferrous ion

= 55.85 x _____ g/L

=

RESULTS & DISCUSSION:

Amount of Ferrous ion present in 1 litre of the given solution = g/L

MAPPING OF PO & PSO: (For all the COs covered by this experiment)

CO No	PO1	PO2
1	2	1

ASSESSMENT:

Particulars	Max. Marks	Marks awarded
Preparation	10	
Conduct of Experiment	30	
Results & Discussion	30	
Viva - Voce	20	
Report	10	
Total	100	
Evaluator's Signature		